

Application Systems

Optics & Optical Coatings

Opto-Mechanics

Bases

Manual Stages

Actuators & Adjusters

Motorized Stages

Light Sources & Laser Safety

Index

Guide

Mirrors

Beamsplitters

Polarizers

Lenses

Multi-Element Optics

Filters

Prisms

Substrates/Windows

Optical Data

Maintenance

Selection Guide

Achromats

Focusing Lenses

f θ Lenses

Objectives

Expanders

Others


Visible spectrum achromats are air spaced achromatic triplets or doublets for lasers in the visible spectrum or white light applications.

The elements are made of crown glass of low dispersion and flint glass of high dispersion.

- These lenses have been optimized for achromatic and spherical aberrations and coma for the 3 wavelengths; blue (486.1nm), green (546.1nm) and red (656.3nm). They are coated with a broadband multi-layer anti-reflection coating for 400 – 700nm.
- Air spaced design allows high power laser applications which includes YAG second harmonic wavelength (532nm).
- The triplets with F-numbers ≥ 2 and doublets with F-numbers ≥ 3 are designed to have each spot size equal to the diffraction limited spot size and very ideal for a Gaussian input beam.




Schematic

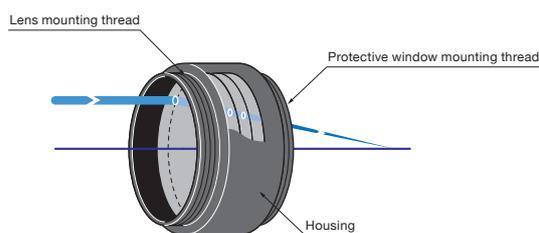

Outline Drawing

(in mm)

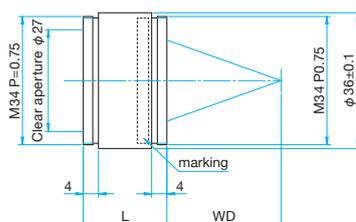

- Tolerance Length $L \pm 0.2$
- Focal length $\pm 2\%$

Specifications

Material	Crown Glass – (Air spaced) – Flint Glass
Material of frame	Aluminum Finishing: Black anodized
Design wavelength	486nm, 532nm, 656nm
Coating	Broadband multi-layer anti-reflection coating
Acceptance angle	$\pm 1^\circ$
Laser Damage Threshold	1J/cm ² (Laser pulse width 10ns, repetition frequency 20Hz)

Guide

- ▶ Please contact our Sales Division for customized achromats. (Customized on size etc.)
- ▶ Protective lens case with rods for mirror holders is available as an option. Please contact us for further information.
- ▶ Please check the “wavelength characteristic of the focal length data” on the Web for the focal lengths of each wavelength.

▶ [WEB Reference](#) [Catalog Code](#) W3078

Attention

- ▶ Since the focal length and working distance of the lens is calculated at 532nm, it will change at other wavelengths due to the refractive index of the material shift.
- ▶ The F number of a lens is calculated by f (effective focal length) / D_e (effective clear aperture). The value represents “Brightness of the lens”. The lower the value, the brighter the lens is.
- ▶ Be sure to wear laser safety goggles when checking optical path and adjusting optical axis.

Specifications

Part Number	Focal length f [mm]	Length L [mm]	Numerical aperture (NA)	Working distance (WD) [mm]
ATL-30-40PY2	40.2	22	0.34	30.1
ATL-30-50PY2	49.4	22	0.27	39.0
ATL-30-60PY2	58.9	22	0.23	49.0
NADL-30-80PY2	80.1	13	0.17	71.8
NADL-30-100PY2	99.8	13	0.14	91.9
NADL-30-150PY2	150.0	12	0.09	142.1
NADL-30-200PY2	199.8	12	0.07	193.1

Compatible Optic Mounts

LHF-M34-30